

## Alkyne Gold Nanoparticles-PEG2K PRODUCT DATA SHEET

# Alkyne Gold Nanoparticles-PEG2K

#### Description

Gold nanoparticles are widely used nanomaterials and generally referred to as colloidal gold in biological research. Colloidal gold markers generally has a particle size between 10 and 100 nm, and will show different colors with the change of particle size. Gold nanoparticles have excellent biocompatibility, rich surface modification properties, and unique optical properties related to the surfactant, shape, size, and structure of the nanoparticles. According to their different characteristics, it can be applied to various fields of biomedicine, such as medical testing, medical imaging, drug delivery, etc.

Cytodiagnostics alkyne functionalized gold nanoparticles are suitable for covalent conjugation of any azide-tagged molecule through Cu(I)-catalyzed click chemistry. These gold nanoparticles allow for rapid generation of gold conjugates for subsequent use in downstream assays such as lateral flow assays and imaging. Our alkyne functionalized gold nanoparticles are available in 11 different sizes ranging from  $5 \sim 100$  nm, are more than 95% spherical and have uniform size distribution (CV < 12%). Products are provided as suspensions at 50 OD.

Abvigen provides a variety of gold nanoparticles, gold nanorods, gold nanocages, gold nanostars, gold nanobipyramids, and other products, the product particle size is optional, the concentration can be customized, the surface can be modified with different groups, and can be appropriately selected according to the customer's use.

For custom sizes, formulations or bulk quantities please contact our customer service department.

Email: info@abvigenus.com

© Abvigen Inc All Rights Reserved

Website: www.abvigen.com Phone: +1 929-202-3014 Email: info@abvigenus.com

#### Characteristics

Composition: Alkyne Gold Nanoparticles-PEG2K

Shape: Spherical

Core diameter: 5 ~ 100 nm

Size dispersity: Coefficient of Variance (CV) < 12%

Polydispersity index (PDI): < 0.200



Size: 1 mL

Amount: OD = 50

Surface functional group: PEG2K-Alkyne

Absorbance (λmax): 510 ~ 570 nm

Alkyne surface density: 1/nm<sup>2</sup>

Buffer: DI Water Form: Suspension

Supplied in H<sub>2</sub>O

**Advantages** 

Monodisperse

Well defined sizes from 5 nm to 100 nm

Precisely engineered surface with an optimized alkyne density

Extensive range of surface functionalities designed for in vitro and in vivo applications

**Application** 

Ideal for development of gold conjugates for use in applications such as blotting, lateral flow assays, light microscopy, and transmission electron microscopy (TEM) among others.

Storage

This product should be stored at 4°C. **DO NOT FREEZE.** If stored as specified, Abvigen Alkyne Gold

Nanoparticles-PEG2K are stable for at least 12 months.

Handling

When stored for a long period of time gold nanoparticles may sediment at the bottom of the vial, which is especially prominent for larger particle sizes. Prior to use, re-suspend the sedimented

particles by thorough mixing until a homogenous solution is obtained.

Note

These products are for R&D use only, not for drug, household, or other uses.




### NPS of Gold Nanoparticles, OD 50

Diameter	Peak SPR	Optical	Wt. conc	Size	Particles/ml	Molarity mol/ml
	Wavelength	density		Dispersity %PDI		
5 nm	515-520 nm	OD 50	2.5 mg/ml	< 20%	1.98E+15	3.28E-09
10 nm	520 nm	OD 50	2.5 mg/ml	< 15%	2.47E+14	4.10E-10
20 nm	524 nm	OD 50	2.5 mg/ml	< 10%	3.09E+13	5.13E-11
30 nm	526 nm	OD 50	2.5 mg/ml	< 6%	9.15E+12	1.52E-11
40 nm	530 nm	OD 50	2.5 mg/ml	< 4%	3.86E+12	6.41E-12
50 nm	535 nm	OD 50	2.5 mg/ml	< 4%	1.98E+12	3.28E-12
60 nm	540 nm	OD 50	2.5 mg/ml	< 4%	1.14E+12	1.90E-12
70 nm	548 nm	OD 50	2.5 mg/ml	< 4%	7.21E+11	1.20E-12
80 nm	553 nm	OD 50	2.5 mg/ml	< 4%	4.83E+11	8.02E-13
90 nm	564 nm	OD 50	2.5 mg/ml	< 4%	3.39E+11	5.63E-13
100 nm	572 nm	OD 50	2.5 mg/ml	< 4%	2.47E+11	4.10E-13

### **Gold Nanoparticles Centrifugation Parameters**

Particle Size	Speed (g)	Time (min)
5 nm	100000	30
10 nm	17000	60 (~ 50% recovery)
20 nm	6500	30
30 nm	4500	30
40 nm	2500	30
50 nm	2000	30
60 nm	1125	30
80 nm	400	30
100 nm	400	30

## TEM of Abvigen gold nanoparticles of different size

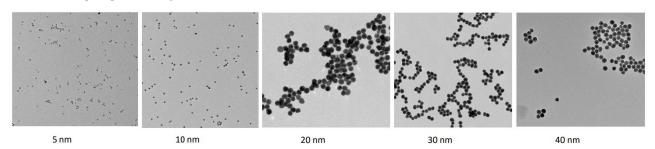

1378 US-206 Ste 6-126, Skillman, NJ USA Tel: 1-816-388- 0112 Fax: 1-888-616-0161 Email: info@abvigenus.com
© Abvigen Inc All Rights Reserved



### **Ordering Information**

Website: www.abvigen.com

Phone: +1 929-202-3014

Email: <a href="mailto:info@abvigenus.com">info@abvigenus.com</a>

1378 US-206 Ste 6-126, Skillman, NJ USA Tel: 1-816-388- 0112 Fax: 1-888-616-0161 Email: info@abvigenus.com © Abvigen Inc All Rights Reserved